

Performer III




reliablecreativesolutions™



Access *and* comfort. All within reach.

Designed to give you flexibility, access and comfort, Performer answers the need for less complexity, lower cost, and exceptional choices for every dollar you invest.

Quiet, smooth, dependable

The chair's thin-profile upholstery allows you to lower the patient to get in close for a comfortable working position. For treating smaller patients and working direct, you can choose Performer's narrow backrest, which gives you exceptional access around the patient's neck and shoulders.

Because patient comfort is just as important as yours, the backrest, seat, and toeboard are designed to cradle the patient and reduce pressure points. For a narrow backrest, there's also the "arm slings" option to accommodate all kinds of body types.



•••
Built to last.
With Performer,
you're assured
comfort and
simplicity.
It's part of what
makes A-dec
A-dec.

•••

Delivery System

Among Performer's distinctive choices is a front- or back-mounting location for the delivery system. The front-mounting location offers more flexibility and provides better ergonomics for those who prefer side delivery and working indirect. The back-mounting location is a good choice if you prefer working direct or need to provide stand-up dentistry.



Assistant's Instrumentation

Performer offers two basic vacuum models: cuspidor with assistant's instrumentation and assistant's instrumentation-only. Both models come standard with a unitized four-position holder, autoclavable saliva ejector, and autoclavable HVE. A syringe can be added as an option.

Performer's extended assistant's arm option, which rotates on the other side of the chair, adds flexibility and provides the solo operator (or assistant) with better access to assistant's instrumentation.







Order Your Samples at a-dec.com/InspireMe




sky blue




indigo



sapphire

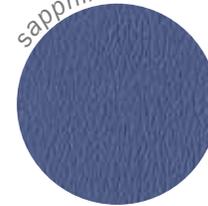

amethyst

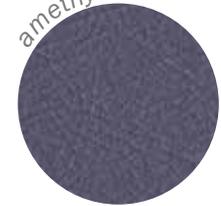

driftwood



black



lemongrass



hazelnut



paprika



tawny



See your authorized A-dec dealer for accurate color samples and the most current product information.



Specifications

Performer III 8000 Dental Chair	Standard	Optional
Thin-profile backrest	●	
Double-articulating headrest	●	
Pivoting armrests		○
Programmable button footswitch	●	
Programmable chair touchpad (on delivery systems)		○
Cast-aluminum baseplate	●	
Hydraulic-lift cantilever system	●	
Accessory power cord		○
Multiple mounting points (for delivery system, cuspidor, light, other options)	●	

Lights

Performer's dental light mounts to a fixed post on the chair and provides both horizontal and vertical adjustment. Choose from the standard Performer halogen 2-axis light or go with the halogen 3-axis A-dec 6300 Light.


Dental Lights	Standard	Optional	Optional
	A-dec Performer Light	A-dec 6300 Light	A-dec LED Light
Rotational axes	Horizontal, vertical	Horizontal, vertical, diagonal	Horizontal, vertical, diagonal
Color temperature	4,500K	5,000K	5,000K
Light pattern	2.4" x 7.1" at 27.6" focal distance (60 mm x 180 mm at 700 mm)	3.3" x 6.3" at 27.6" focal distance (85 mm x 160 mm at 700 mm)	3.8"x 5.7" at 27.6" focal distance (95 mm x 145 mm at 700 mm)
Light intensity	Composite: 8,000 lux (743 fc) High: 11,500 lux (1068 fc)	Composite: 8,000 lux (743 fc) Medium: 20,000 lux (1858 fc) High: 24,000 lux (2230 fc)	Low: 15,000 lux (1394 fc) Medium: 25,000 lux (2323 fc) High: 30,000 lux (2787 fc) Cure-safe Illumination: 25,000 lux (2323 fc) yellow light
BTUs per hour	180	325	24

Chair-mounted Accessory Tray

Performer's unique accessory tray mounts on either side of the chair and rotates for chair-side or trans-thorax positioning. Capable of supporting up to 12 lbs (15.4 kg), this convenient tray can hold instruments, serve as a periodontal charting surface, or support accessories such as a tabletop scaler.

Legendary A-dec reliability

You can count on A-dec equipment to work when you do. Reliable components on the inside and attention to detail everywhere else. With A-dec, you get equipment that looks beautiful and works consistently. Day after day. Year after year.

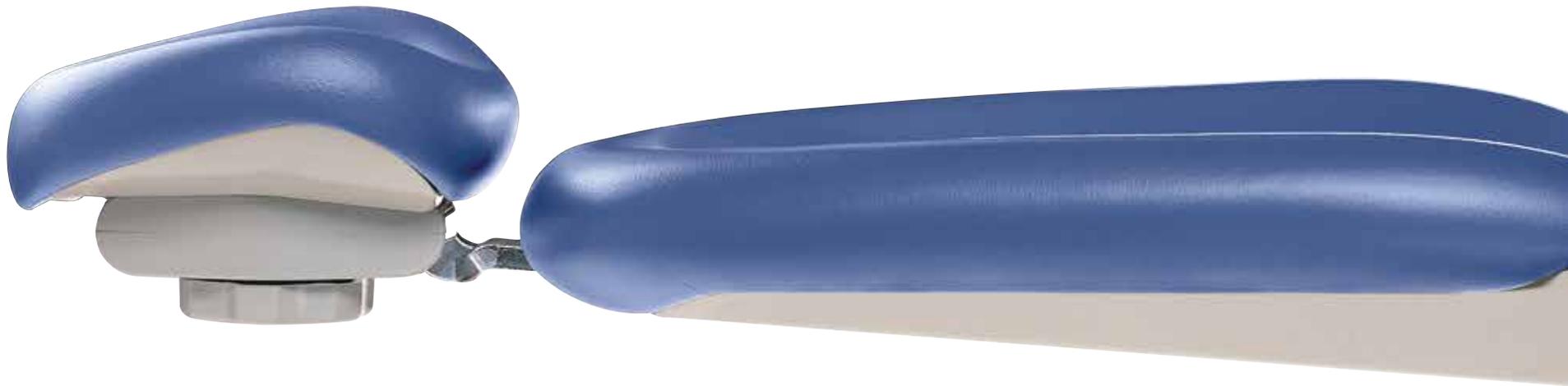

A-dec Headquarters

2601 Crestview Drive
Newberg, Oregon 97132 USA
Tel: 1.800.547.1883 within USA/CAN
Tel: 1.503.538.7478 outside USA/CAN
Fax: 1.503.538.0276
a-dec.com

A-dec Australia

Tel: 1.800.225.010 within AUS
Tel: +61.(0)2.8332.4000 outside AUS

A-dec China

Tel: +1.503.538.7478

A-dec United Kingdom

Tel: 0800.ADEC.UK (2332.85) within UK
Tel: +44.(0).24.7635.0901 outside UK

A-dec, the A-dec logo, A-dec 500, A-dec 300, Cascade, Cascade Master Series, Century Plus, Continental, Decade, ICX, ICV, Performer, Preference, Preference Collection, Preference ICC, and Radius are trademarks of A-dec Inc. and are registered in the United States and other countries. A-dec 400, A-dec 200, Preference Slimline, and reliablecreativesolutions are also trademarks of A-dec Inc. None of the trademarks or trade names in this document may be reproduced, copied, or manipulated in any manner without the express, written approval of the trademark owner.

Certain touchpad symbols are proprietary to A-dec Inc. Any use of these symbols, in whole or in part, without the express written consent of A-dec Inc., is strictly prohibited.

The colors and products shown may vary from the actual materials and are subject to change without notice. Please see your authorized A-dec dealer for accurate color samples and the most current product information.