



# Clesta II

**Clesta II treatment centres represent outstanding flexibility** and the very best of Belmont's hygiene-conscious design. Air and Electric options across the comprehensive range of delivery systems mean that any Practice or Hospital can fulfill its exact requirements. The improved comfort of Clesta II's chair is an ideal base on which to build your dental system.

## CLESTA II-CM<sup>1</sup>

with unit-mounted 720S operating light<sup>2</sup>

Holder type **Air** or **Electric** versions

Continental Rod type **Electric** version



## CLESTA II 'Flexible' Treatment Centre (a) or (b)<sup>3</sup>



**CLESTA II-MC<sup>4</sup>**  
Air or Electric<sup>5</sup> versions



**CLESTA II-M<sup>4</sup>**  
Air or Electric<sup>5</sup> versions

<sup>1</sup> Available with left- or right-handed factory installation / <sup>2</sup> 'Alternative' Clesta II-CM available without unit-mounted 720S operating light  
<sup>3</sup> Without Assistant's console / <sup>4</sup> Available separately / <sup>5</sup> 720S operating light, see page 23 / Clesta II operating light options, see page 23

# Clesta II over-patient delivery :

*Versatile, balanced and precise* delivery systems providing total working comfort between the 9 and 1 o'clock positions.

## CLESTA II-CM HOLDER TYPE

Rotating spittoon (180°) with cup-filler and bowl-flush

Smooth, balanced arm movement with pneumatic brake

### 'E' version key features ~

- LED digital display
- Preset micro motor
- Torque control with Auto-Reverse
- First Priority handpiece selection

Height-adjustable, easily positioned and zero drift for **comfortable operation** at any position between 9 and 1 o'clock

Removable rubber mat and extensive console table

Ultra-slim backrest and raised sides **improve Operator's working comfort**



# 'Holder' or 'Continental Rod' :

# Clesta II

*Ergonomic and hygienic* features are standard on all Clesta II treatment centres.

## CLEST II-CM CONTINENTAL ROD TYPE

Easy-clean, detachable parts including removable spittoon bowl

Seamless upholstery

Seamless, easy-clean upholstery offers optimum hygiene levels

Control chair movements and spittoon functions from extendable, rotating Assistant's console

Compact control panel with LED display

Soft rubber, non-slip instrument mat is removable to aid hygiene maintenance

*Cuspidor unit is positioned lower than usual along the seat to increase access for an Assistant at the 3 o'clock position*

# Ergonomics & comfort :


*Comfort class for Dentist and Patient.*

Over-patient versatility for comfortable single- and four-handed working positions



## **Chair movements**

*are synchronised and smooth, eliminating the anxiety and insecurity that is common among Patients*

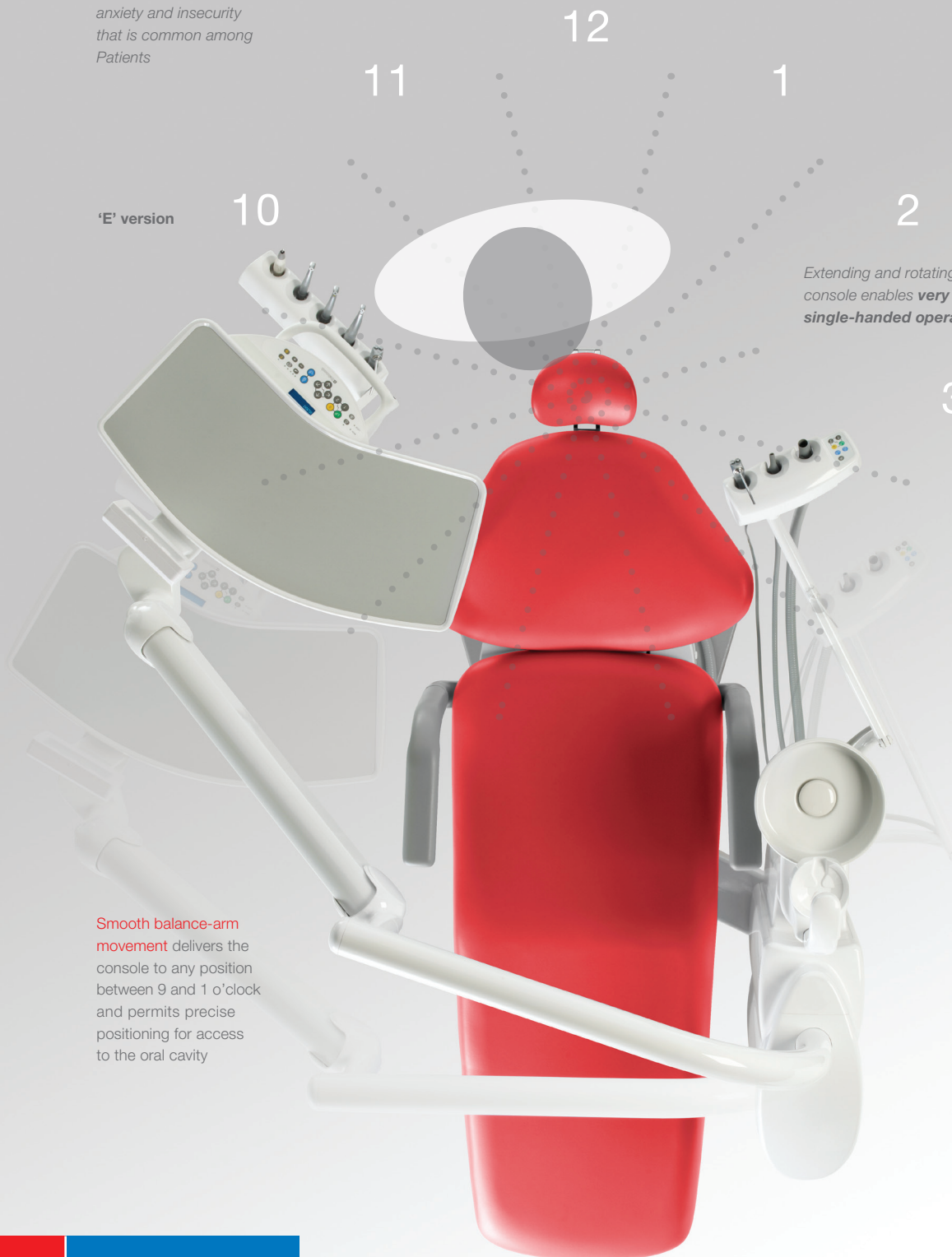

'E' version

*Extending and rotating the Assistant's console enables **very comfortable single-handed operation***

**Smooth balance-arm movement** delivers the console to any position between 9 and 1 o'clock and permits precise positioning for access to the oral cavity

# cuspidor & control :

# Clesta II

*Excellent Assistant access at the 3 o'clock position and easy-clean, hygiene-conscious design.*

**SP-A** Easy-clean, detachable cup-filler holder, spittoon bowl and vacuum tubes

Rotating spittoon (180°)



Spittoon bowl options:  
blue/green/pink glass  
(porcelain as standard)

Assistant's console with chair, spittoon and light<sup>1</sup> controls



'E' version control panel includes instrument, chair, spittoon and light controls



'A' version control panel includes chair and spittoon controls

**Control all chair movements** by hand or foot: by hand from both the Operator's and Assistant's console; by foot using the foot control. Spittoon functions and ON/OFF light control<sup>1</sup> are also activated from both the Operator's and Assistant's console.

<sup>1</sup> Light controls on 'E' version only

# Clesta II 'Flexible' treatment centres :

*Mobile Cart and Cabinet Module delivery systems offer space-saving qualities and the opportunity to configure the Practice's treatment centre exactly as desired.*

**CLESTA II 'FLEXIBLE' TREATMENT CENTRE**

**CLESTA II-MC – 'A' or 'E' version**

Height-adjustable Operator's console

Extendable, rotating Assistant's console

Removable rubber mat



'A' version Mobile Cart

Instruments foot control

Smooth-running castors

equipped & installed :

**Clesta II**

*Clesta II treatment centres can be tailored exactly to suit a busy Practice's requirements.*

**700 SERIES LIGHT** – Unit- or ceiling-mounted<sup>1</sup>

Retractable Patient mirror



ON/OFF Touchless sensor and soft-push button (with light intensity indicator) for Composite Curing Safe Mode and two light intensity settings: Low – 18,000 Lux and High – 28,000 Lux. **4200 Kelvin**

**CLESTA II-MC** and **CLESTA II-M** 'A' and 'E' versions



**SP-A** cuspidor unit with extendable and rotating Assistant's console



**SP-B** cuspidor unit without Assistant's console



Assistant's console with chair, cuspidor and light<sup>2</sup> controls



Continental Rod type Operator's console ('E' version)

<sup>1</sup> NDL-PG unit-mounted operating light also available / <sup>2</sup> Applies to 'E' versions only and when using a 720S Operating light  
Ceiling-mounted NDL and 700 Series lights also available, see 'Operating Lights' on pages 48-51